


USER GUIDE


FOR THE SCANNING AND ANALYSIS SOFTWARE

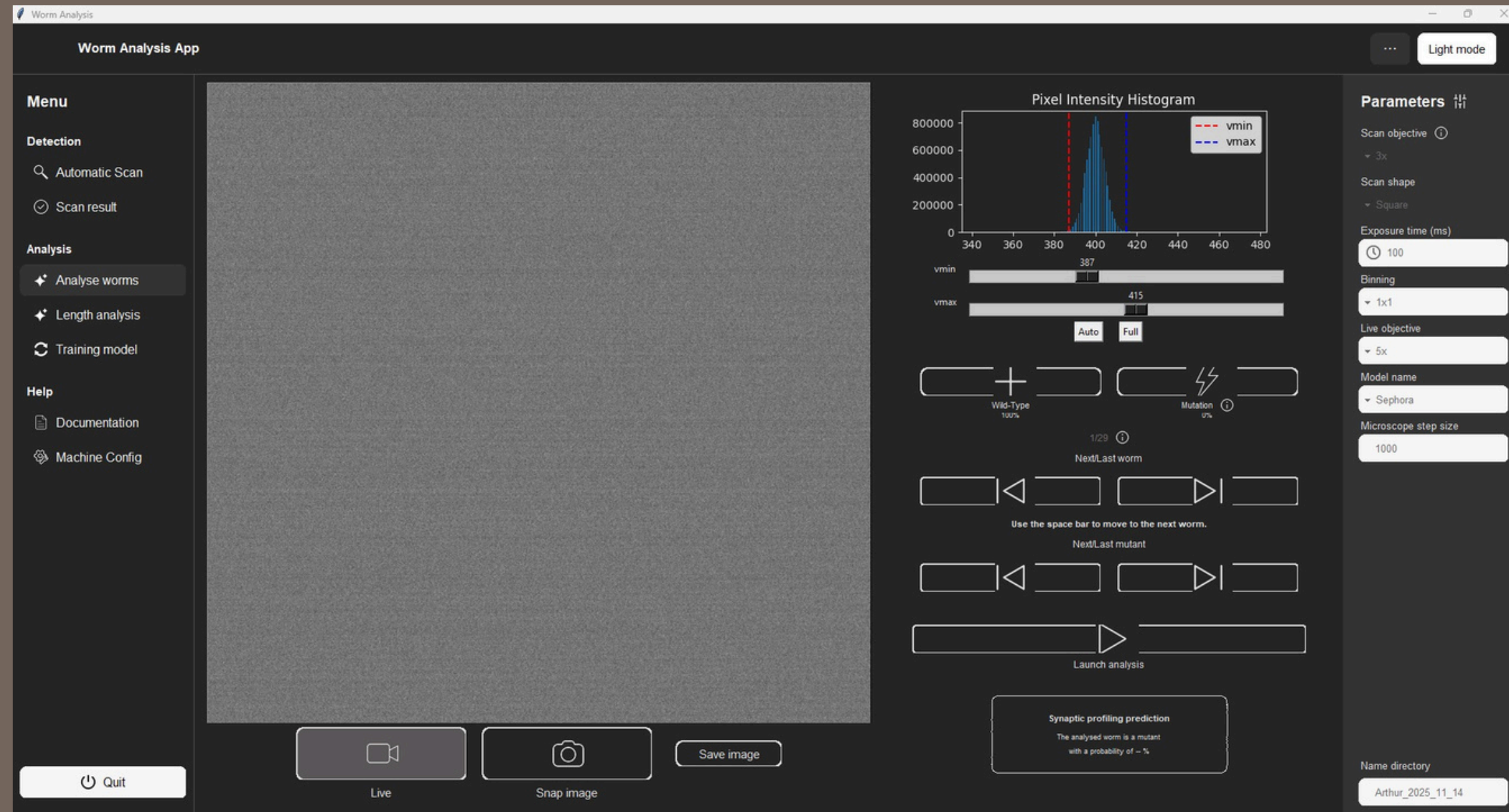

ARTHUR SKOWRONEK

arthur.skowronek13@gmail.com

Open and close the settings bar

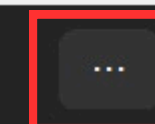

Light mode

Menu

Detection

Automatic Scan

Scan result

Analysis

Analyse worms

Length analysis

Training model

Help


Documentation

Machine Config

Quit



Menu section :
Click on the
desired title to
navigate to a
different
page.

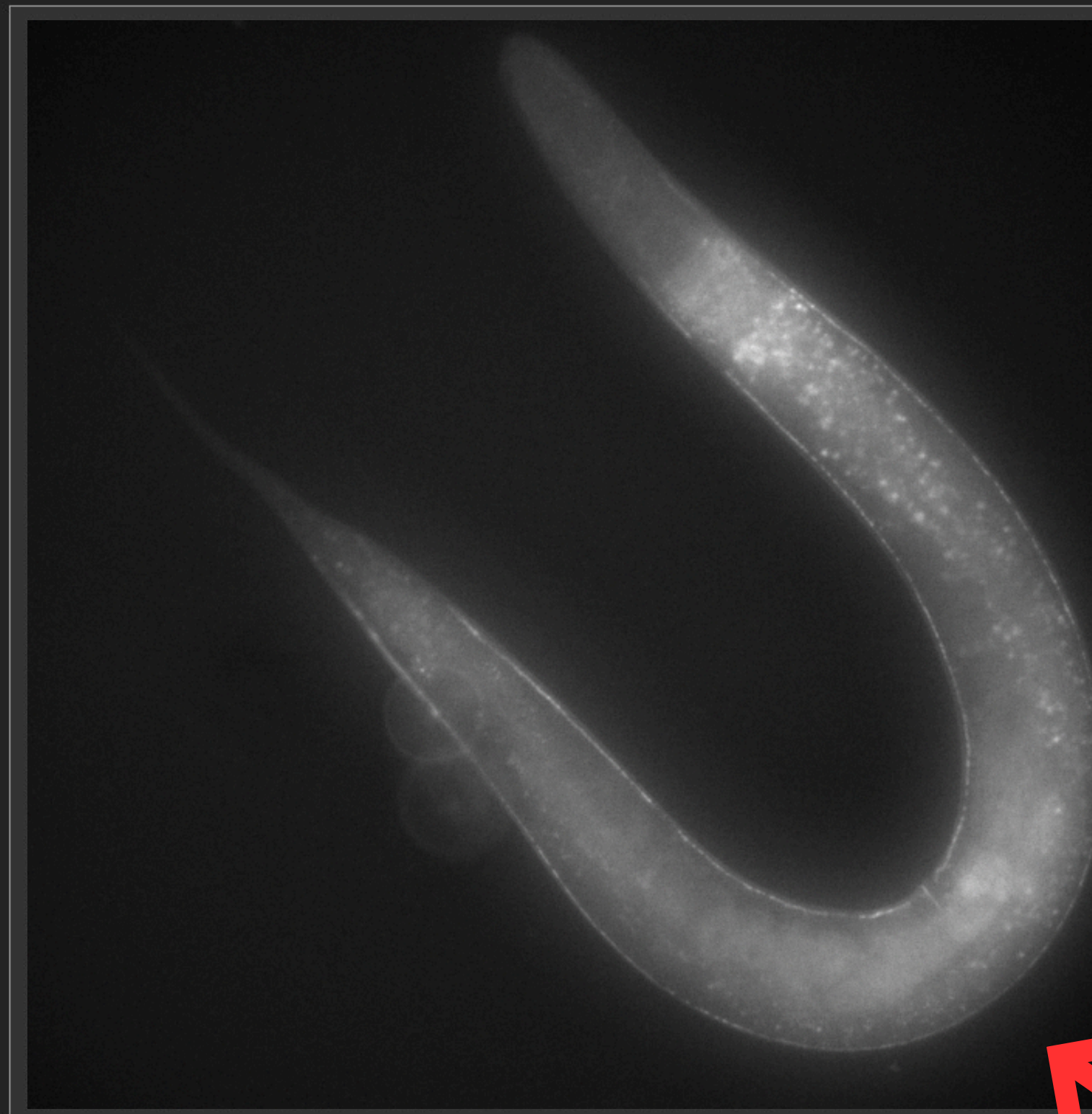

Live image

Start the
scanning
program



Launch scan

Parameters

Scan objective

3x

Scan shape

Square

Exposure time (ms)

100

Binning

1x1

Live objective

5x

Model name

Sephora

Microscope step size

1000

Name directory

Arthur_2025_11_14

How to scan ?

Worm Analysis App

Light mode

Menu

Detection

- Automatic Scan
- Scan result

Analysis

- Analyse worms
- Length analysis
- Training model

Help

- Documentation
- Machine Config

Parameters

Scan objective ⓘ

- 3x
- 2x
- 3x
- 4x
- 5x
- 6x
- 7x

100

Binning

- 1x1

Live objective

- 5x

Model name

- Sephora

Microscope step size

1000

Name directory

Arthur_2025_11_14

1. Use the right objective :

- Macrozoom : x2 with the dial 3
- Nikon : x4

2. Select the used objective (corresponds to the 3 dial for the Macrozoom)

3. Focus the image

4. Move the microscope to the lower right corner

5. Launch the scan

Quit

Launch scan ⓘ

NOTES


1. If you are looking at tiny worms : you can use a higher magnification to do the scan.

For the Macrozoom : ALWAYS use the x2 objective.
But you can change the dial.

For the Nikon : Use whichever objective you prefer
In both cases : Don't forget to update the value in the parameters panel in "Scan Objective"

2. You can set the scan shape to either square or rectangle.

Scan shape : 

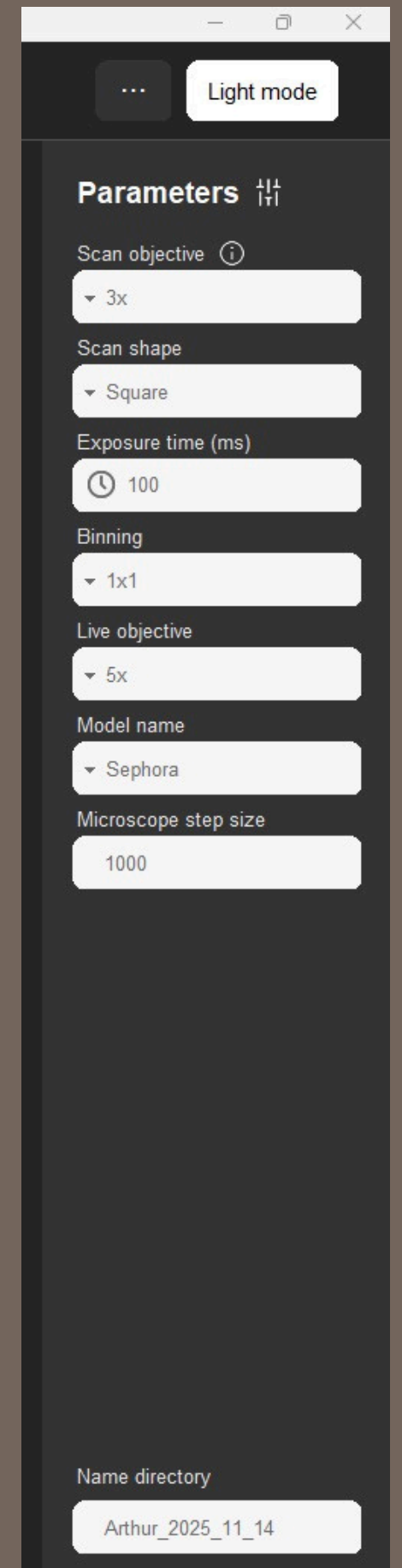

Parameters

Scan objective ⓘ

▼ 3x

Scan shape

▼ Square

Exposure time (ms)

🕒 100

Binning

▼ 1x1

Live objective

▼ 5x

Model name

▼ Sephora

Microscope step size

1000

Name directory

Arthur_2025_11_14

Worm Analysis App

...

Light mode

Menu

Detection

Analysis

Help

Automatic Scan

Scan result

Analyse worms

Length analysis


Training model

Documentation

Machine Config

Quit

After performing the scan : clean the results by adding undetected worms or removing errors

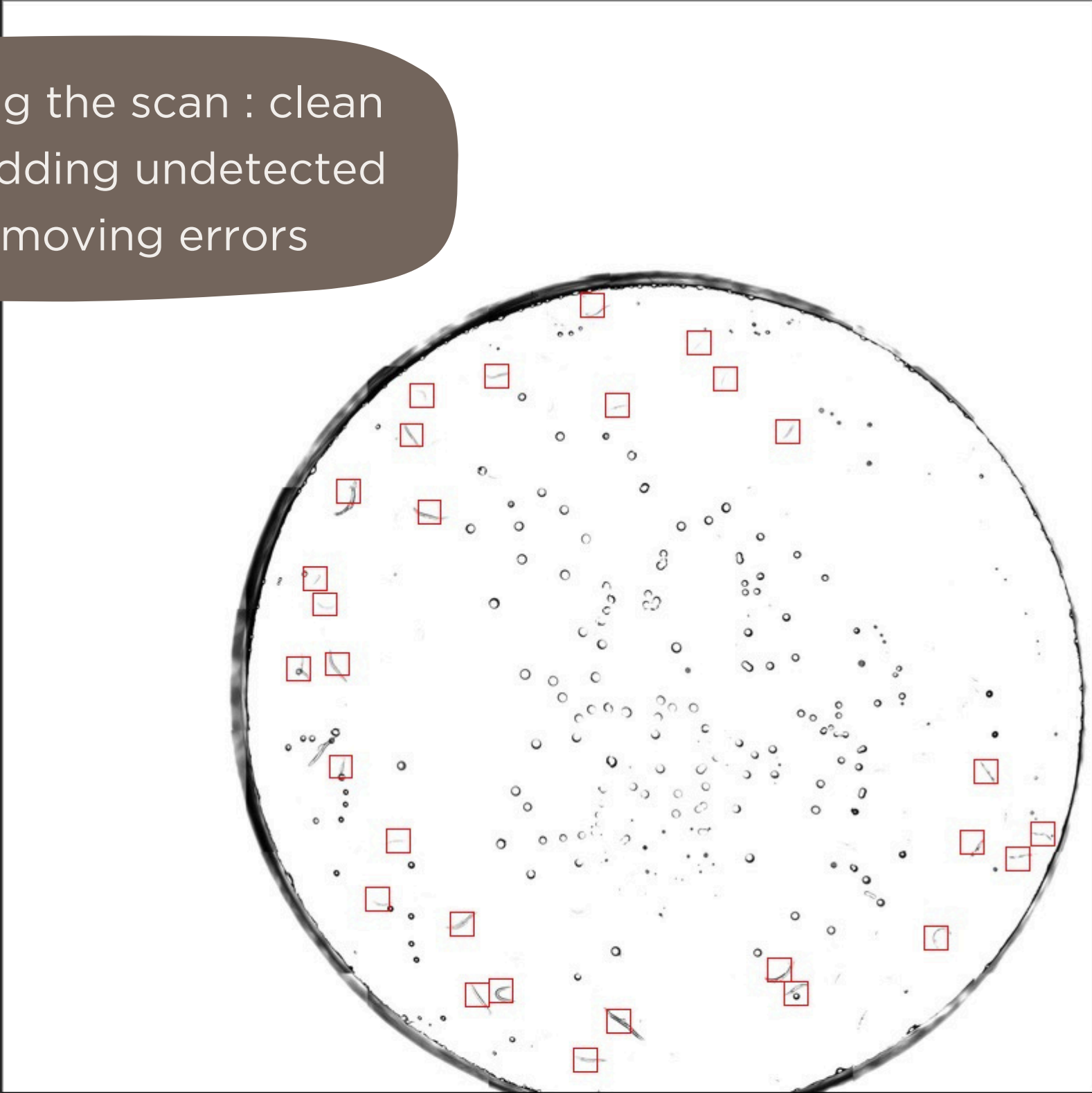

+

+

×

- Add worm** : you only have to click on the image where you want to add a worm
- Remove worm** : click on the worm you want to delete. Or drag your mouse on the worms
- Move worm** : click on a red rectangle and drag it to the desired position

Select the mode you want to use to clean the worm detection

Parameters

Scan objective

3x

Scan shape

Square

Exposure time (ms)

100

Binning

1x1

Live objective

5x

Model name

Sephora

Microscope step size

1000

Name directory

Arthur_2025_11_14



Menu

Detection

- Automatic Scan
- Scan result

Analysis

- Analyse worms
- Length analysis
- Training model

Help

- Documentation
- Machine Config

Parameters

- Scan objective 3x
- Scan shape Square
- Exposure time (ms) 100
- Binning 1x1
- Live objective 5x
- Model name Sephora
- Microscope step size 1000

Counter : Displays the ID of the worm being viewed over the total number of detected worms

1/29

Next/Last worm



Use the space bar to move to the next worm.

Tell the microscope to move to the next (or last) worm position

You can save an image regardless of the mode used (live or snap). The image will be saved in the folder defined here: (the folder will be created automatically if it does not exist)

Name directory

Arthur_2025_11_14



Menu

Detection

- Automatic Scan
- Scan result

Analysis

- Analyse worms
- Length analysis
- Training model

Help

- Documentation
- Machine Config


Parameters

- Scan objective 3x
- Scan shape Square
- Exposure time (ms) 100
- Binning 1x1
- Live objective 5x
- Model name Sephora
- Microscope step size 1000

You can tag the worm being viewed as having a mutation or not



By hovering the info icon you can see the list of the ID of the worm being tag as Mutant

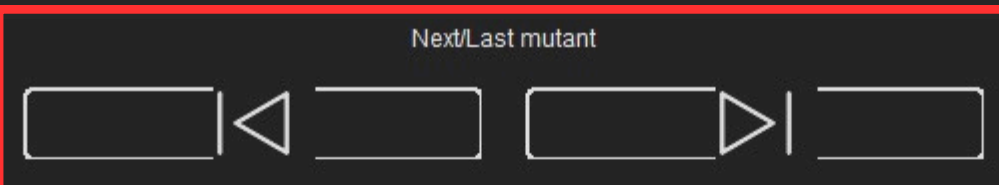

Tell the microscope to move to the next (or last) worm with a mutation

Worm Analysis App

...

Light mode

Menu



Detection

 Automatic Scan Scan result

Analysis

 Analyse worms Length analysis Training model

Help

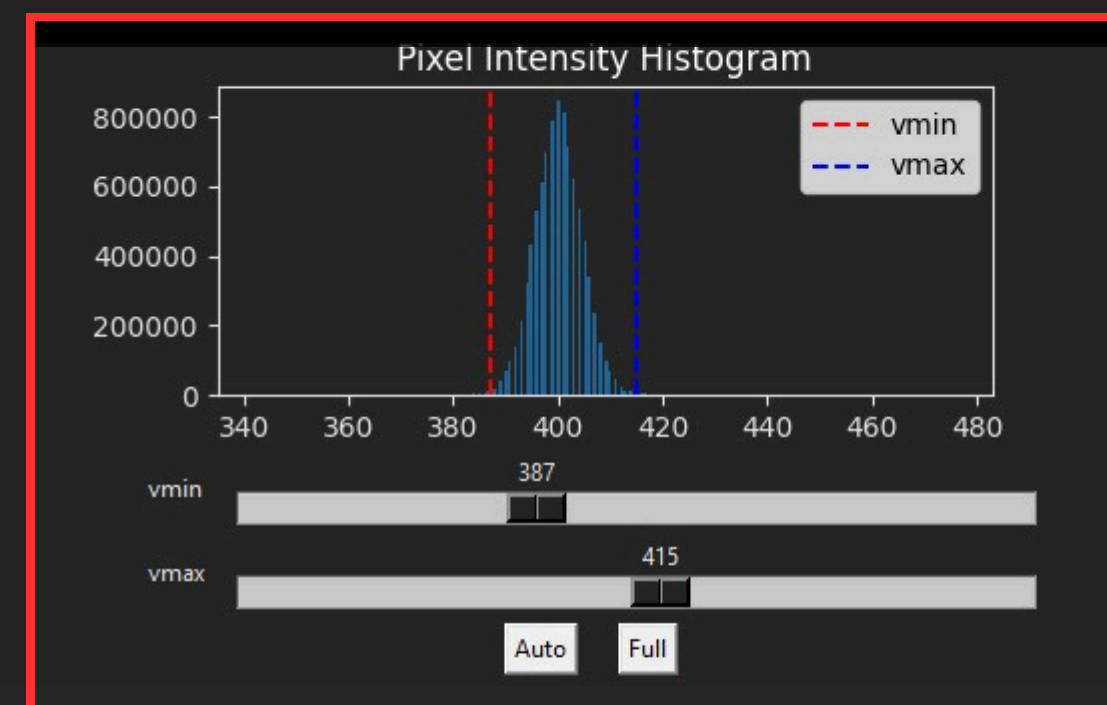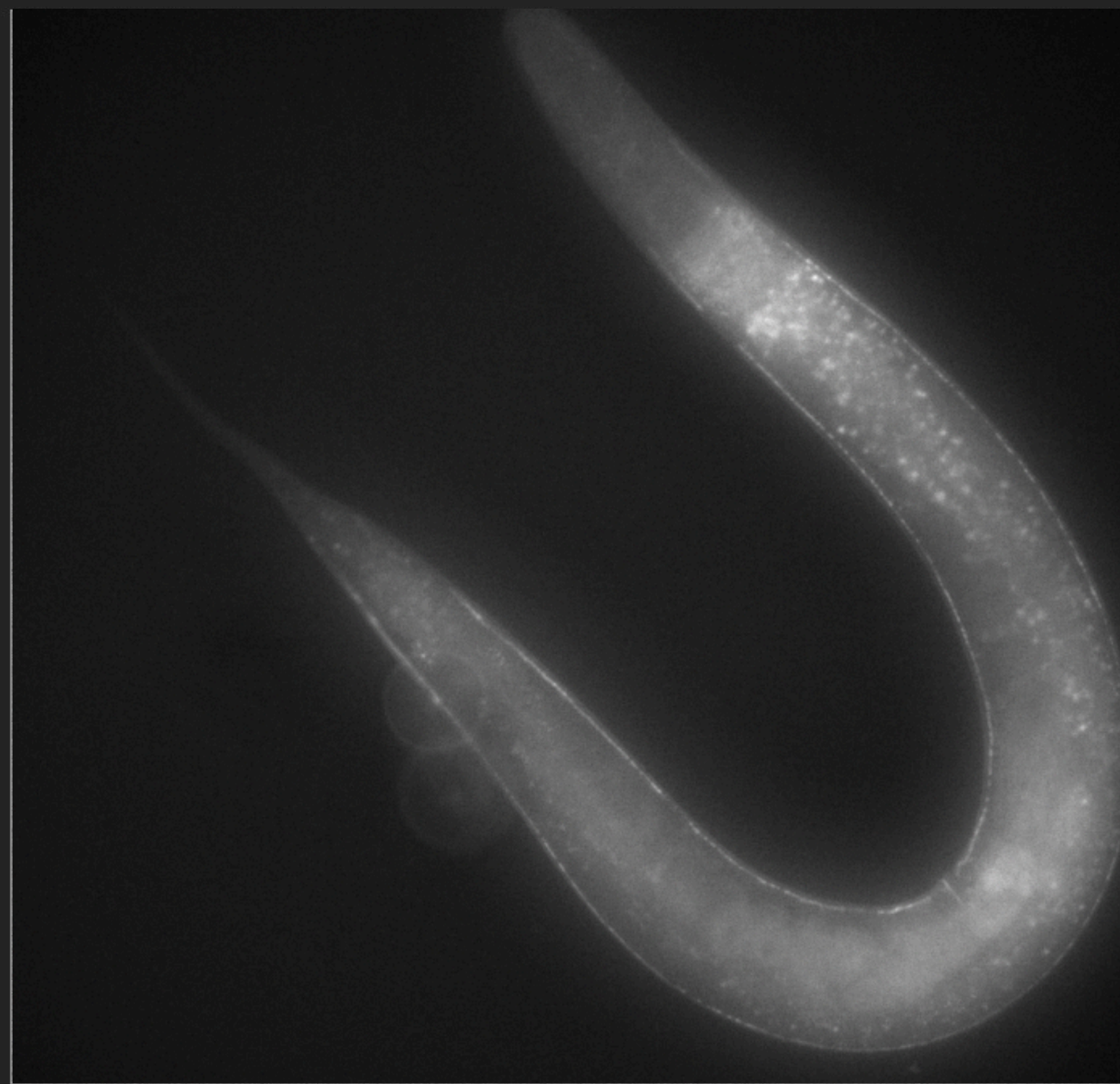 Documentation Machine Config

Change the image displayed
with the sliders and buttons

Parameters Scan objective 

▼ 3x

Scan shape

▼ Square

Exposure time (ms)

 100

Binning

▼ 1x1

Live objective

▼ 5x

Model name

▼ Sephora

Microscope step size

1000

Name directory

Arthur_2025_11_14

 Quit

Live



Snap image

Save image

Worm Analysis App

...

Light mode

Menu

Detection

 Automatic Scan Scan result

Analysis

 Analyse worms Length analysis Training model

Help

 Documentation Machine ConfigParameters Scan objective 

▼ 3x

Scan shape

▼ Square

Exposure time (ms)

 100

Binning

▼ 1x1

Live objective

▼ 5x

Model name

▼ Sephora

Microscope step size

1000

Select the wanted AI model :

Launch the analysis to get a
prediction about the worm mutation

Launch analysis

Synaptic profiling prediction

The analysed worm is a mutant
with a probability of – %

Name directory

Arthur_2025_11_14

 Quit

Live



Snap image

Save image

Worm Analysis App

...

Light mode

Menu

Detection

- Automatic Scan
- Scan result

Analysis

- Analyse worms
- Length analysis
- Training model

Help

- Documentation
- Machine Config

Live image : If you click on it, the microscope move to center to the clicked point

Quit

Live

Snap image

Save image

On the Nikon microscope, select the used objective during the analysis (necessary if you want to click on the live image to move the microscope)

You can use the arrow keys on your keyboard to move the microscope. The size of a step is defined here :

Parameters

Scan objective

- 3x

Scan shape

- Square

Exposure time (ms)

100

Binning

- 1x1

Live objective

- 5x

Model name

- Sephora

Microscope step size

1000

Name directory

Arthur_2025_11_14

Worm Analysis App

...

Light mode

Menu


Detection

 Automatic Scan Scan result

Analysis

 Analyse worms Length analysis Training model

Help

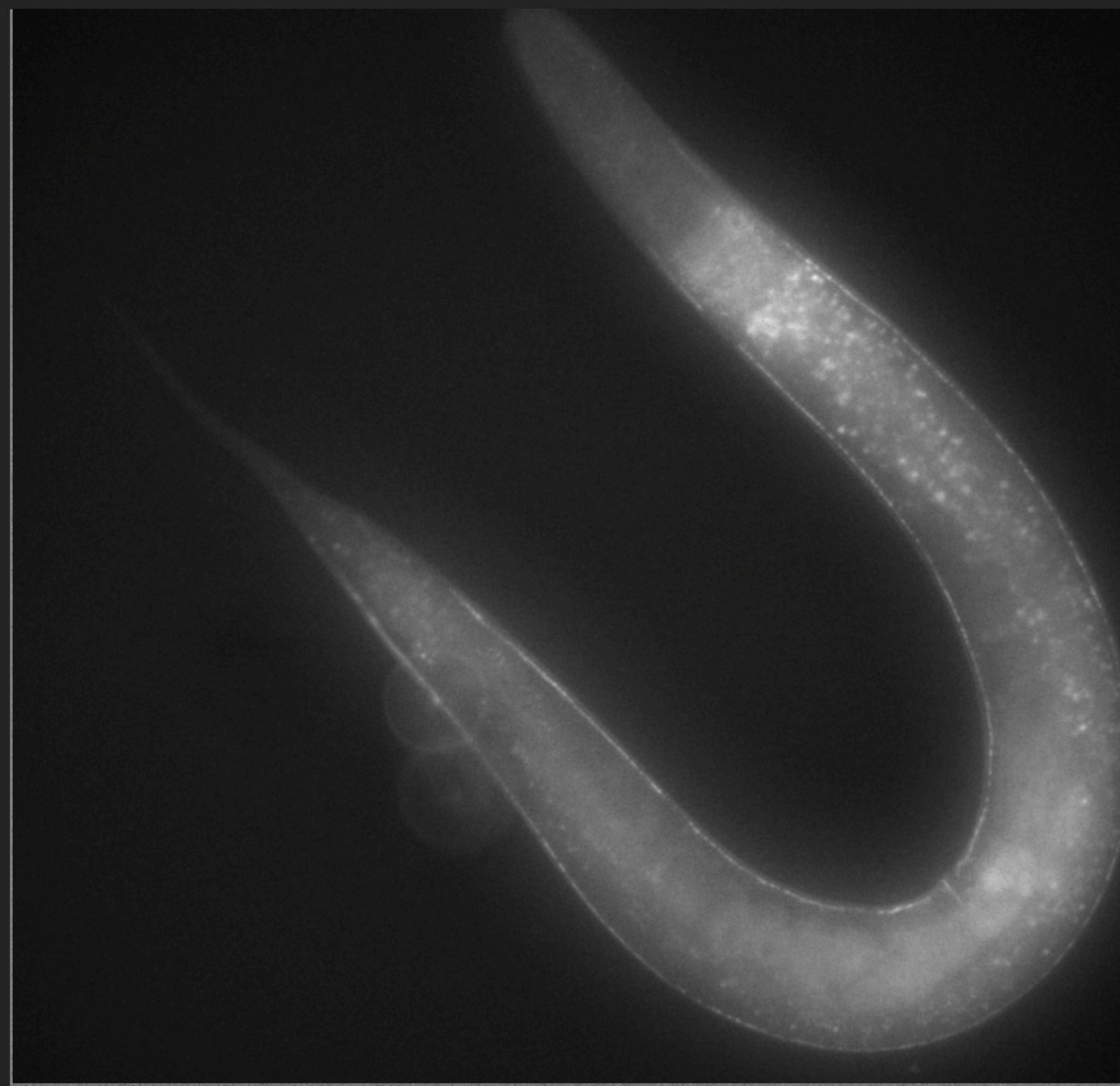 Documentation Machine Config Quit


Live



Snap image

Save image

Pixel Intensity Histogram

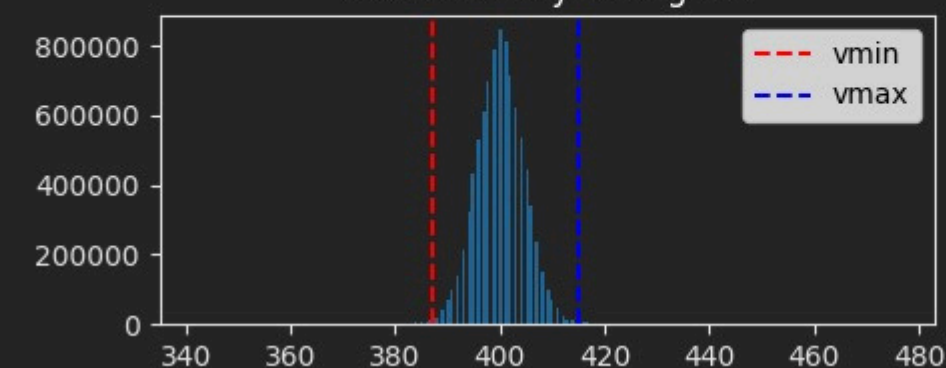

vmin

vmax

Auto

Full

Wild-Type 100%

Mutation 0% 

1/29 

Next/Last worm



Use the space bar to move to the next worm.

Next/Last mutant



Launch analysis

Synaptic profiling prediction

The analysed worm is a mutant
with a probability of – %

Parameters Scan objective 

3x

Scan shape

Square

Exposure time (ms)

100

Binning

1x1

Live objective

5x

Model name

Sephora

Microscope step size

1000

Name directory

Arthur_2025_11_14

Worm Analysis

Worm Analysis App

...

Light mode

Menu

Detection

Analysis

Help

Quit

Automatic Scan

Scan result

Analyse worms

Length analysis

Training model

Documentation

Machine Config

1. Give a name to your model
(either create a new one or
improve an existing one)

2.Drop images in the
following folders to
train a new model

3. Launch the training

Training a model

Step 1 : Choose or create a model

Step 2 : Add images

Step 3 : Train the model

How to train a new AI model
for the mutation detection ?

4.To use this new
model, you can select
it here

Parameters

Scan objective

3x

Scan shape

Square

Exposure time (ms)

100

Binning

1x1

Live objective

5x

Model name

Sephora

Microscope step size

1000

Name directory

Arthur_2025_11_14

WARNING : a AI model can't be moved
between computers. You have to create
it on the computer you will use

You don't necessarily need to have
the camera turn on to train a model

Worm Analysis App



Light mode

Menu

Detection

Automatic Scan

Scan result

Analysis

Analyse worms


Length analysis

Training model

Help

Documentation

Machine Config

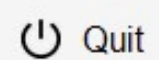

Quit

Parameters

Scan objective

▼ 3x

Scan shape

▼ Square

Exposure time (ms)

100

Binning

▼ 1x1

Live objective

▼ 5x

Model name

▼ Sephora

Microscope step size

1000

Name directory

Arthur_2025_11_14

Be careful when you change these values.
Once you have configured them for your system's operation, do not change them again.
Incorrectly configured values will prevent the program from working properly.

Change the scan length

26000

45000

Manage the scan area (default values : 26 000 ; 45 000)

2

3

4


5

6

7

Manage the magnifications on your microscope

List all the objective size available on the microscope

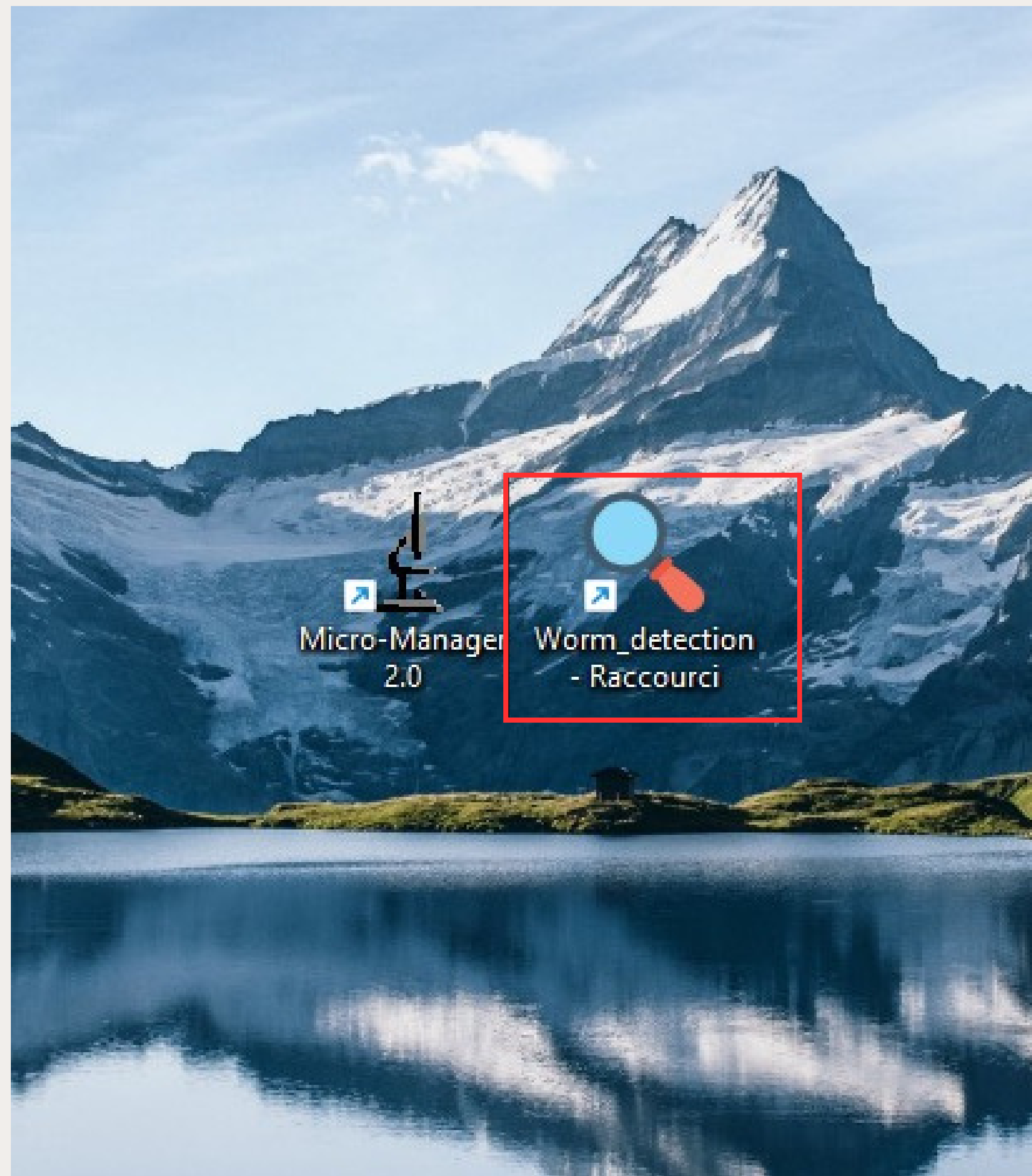

Click on this app to open the software.
The launching can take 10-20s.